

Installation cube


Weak ceiling structure? Expensive installation?

With the Arcoma Intuition Installation cube it is now possible to install a ceiling suspended x-ray tube support in rooms with weak ceiling structure or even in rooms without a ceiling! The Installation cube will give you the full flexibility of a ceiling suspended x-ray tube support where you until now have been limited to other solutions.

Also consider the Installation cube as a tool for a cost effective installation. Add the floor mounted Installation cube instead of doing a costly reconstruction of the ceiling. You will have full control of cost and time of the installation – no surprises!



Installation cube

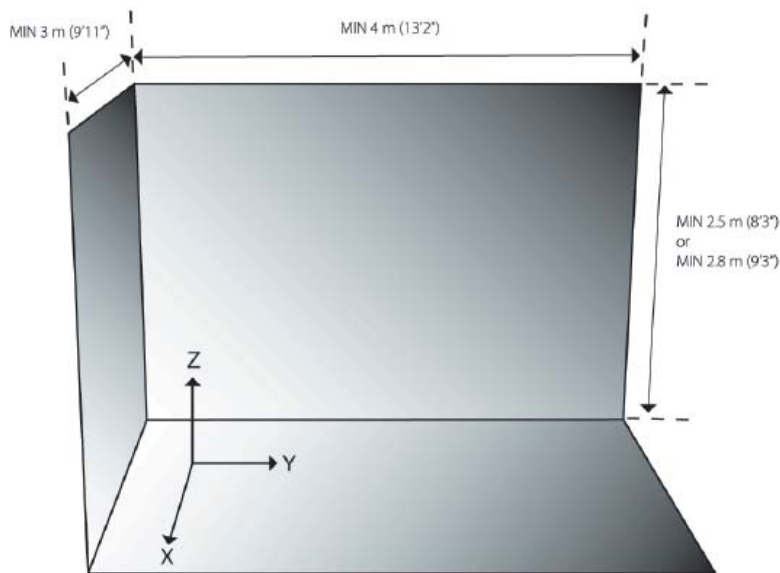

Two different versions

The Installation cube is available in two different versions. Version 1 is intended for really small rooms requiring only 2.5 m (8'3") ceiling height. Version 2 gives you an increased stroke in the z-direction for rooms with ceiling height of 2.8 m (9'3") or more.

Usable stroke

The length of the stroke of the ceiling suspended x-ray tube:

- X: 1850 mm (6'1")
- Y: 2750 mm (9'1")
- Z: 1450 mm (4'9"), Version 1: Ceiling height 2,5m (8'3")
1700 mm (5'7"), Version 2: Ceiling height 2,8m (9'3")

Focuspoint to floor: 1

approximately 300 – 1750 mm (1' – 5'9"), Version 1: Ceiling height 2,5m (8'3")
approximately 300 – 2000 mm (1' – 6'7"), Version 2: Ceiling height 2,8m (9'3")

1) Depending on the selected tube